

Vapocid[®]-Tray System



A convincingly simple and cheap solution for ideal sterile instrument storage. Triple benefit – the system can be used in three different ways: as a working tray, a washing tray and as sterile packaging.




There are many cases where it can be technically difficult or costly to meet modern hygiene requirements – but this is where the patented Vapocid®-Tray offers a real solution, because it can be used in three ways: as a washing tray (e.g. in thermo-disinfectors), as a normal working tray, and as sterile packaging within a tray system. Thanks to these features, the highest economic and ecological requirements can be attained.

(See report in "Schweiz. Monatszeitschr. Zahnmed.[Swiss Monthly Dental Journal]", pages 1061 to 1072, vol. 109, 10/99).



Individual assembly means that each working tray can be optimally equipped. The plastic fixtures will also hold the instruments in a diagonal or slanted position. A washing cover can be used for washing and disinfecting, eliminating the need for an extra washing basket. The Vapocid®-Tray can be stacked so as to save space during transport and storage.

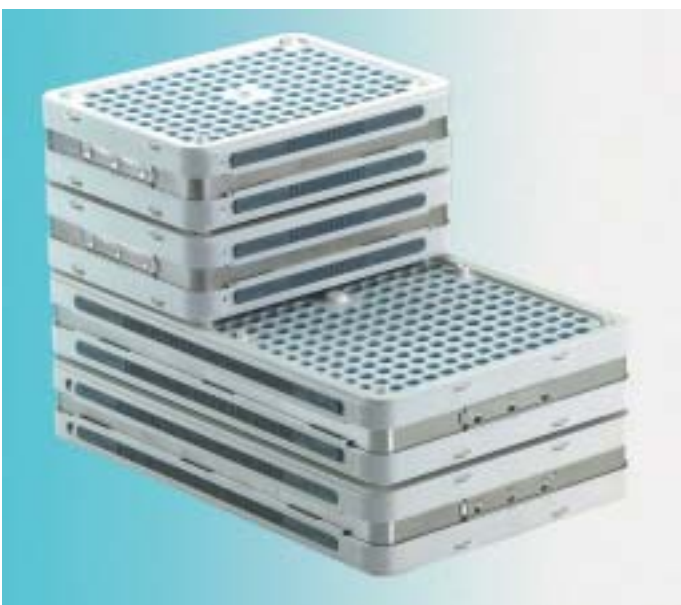

The Vapocid®- Tray conforms to CE and it can be supplied in a variety of sizes and materials.


Part No. 38501 -1/1- Standard tray (288 x 189 x 47 mm)

*Part No. 38502 -1/2- Standard tray (189 x 144 x 47 mm)
made entirely of special steel*

Part No. 38525 -1/1- Standard tray, made of special steel

*Part No. 38523 -1/2- Standard tray, made of special steel
both with plastic covers*

The plastic covers (top + bottom) can always be replaced with special stainless steel parts at a later stage.

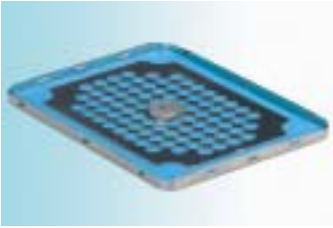


Structure of the Vapocid®-Tray special steel covers (top + bottom) with filter and filter holder



Special steel working, washing, disinfecting and sterile tray with two label holders for inscriptions



Special steel sealing clips with labelling area for batch number and indication label

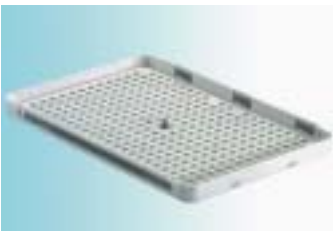

*Structure of covers (top + bottom) to maintain sterility
Plastic cover with sealing groove, slide lock and screw for filter holder. This cover is also used for washing.*



Filter with welded-on silicone seal. Suitable for multiple use. Waste disposal costs are cut substantially because disposable packaging is no longer needed (thus environmental pollution is also reduced).



Special steel filter holder to protect the filter (speeds up the sterilisation procedure).



Plastic instrument holder (easy to fit)



Seal for one-time use (example) on request only




Special steel label for for inscriptions.



Special CNC machine to apply the sealing compound



Manufacturing filters for the patented Vapocid®-Tray System

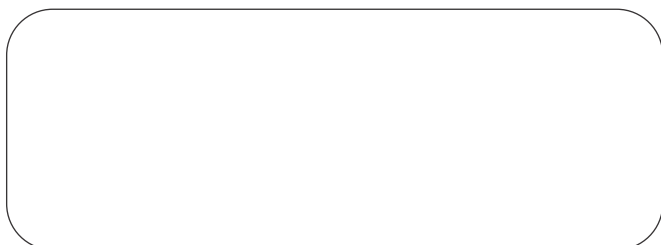

OMB MEDITECH AG
Industrie Hofacker, CH-9606 Bütschwil · Switzerland
Tel: +41-71-982-7201
Fax: +41-71-982-7202
URL: www.omb-meditech.com
eMail: sales@omb-meditech.com